


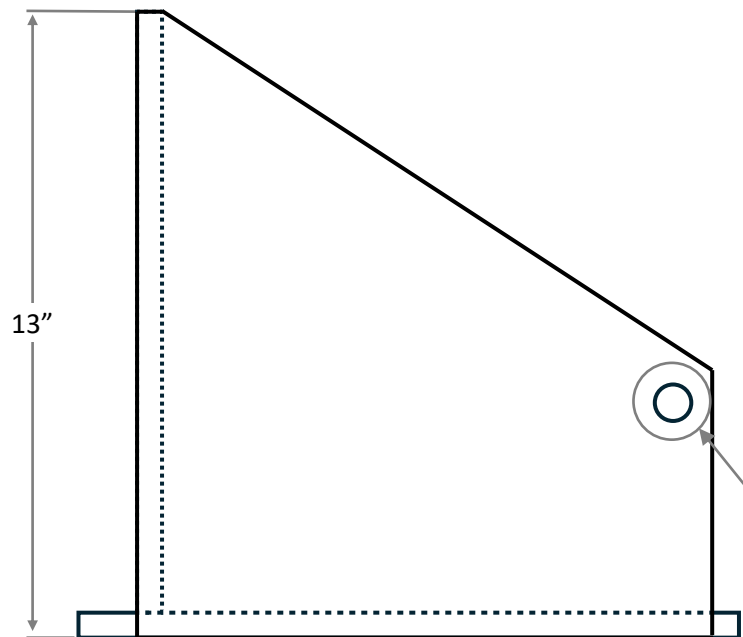
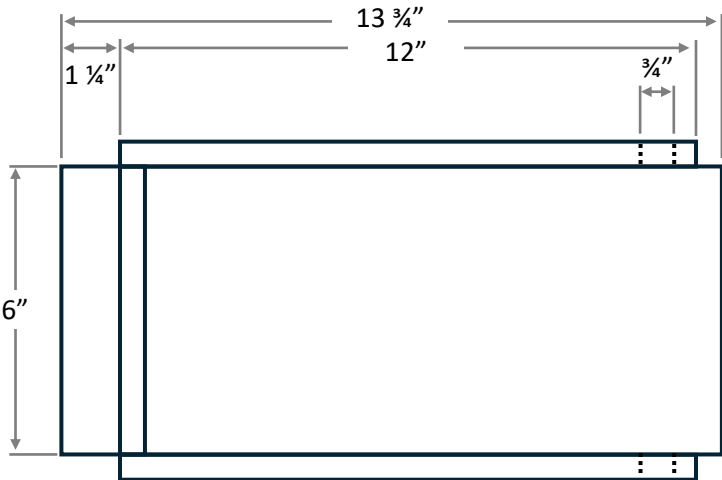


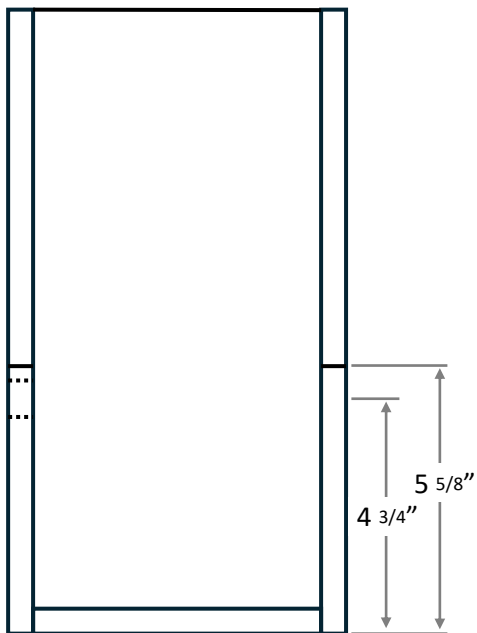


See detail for
axis assembly

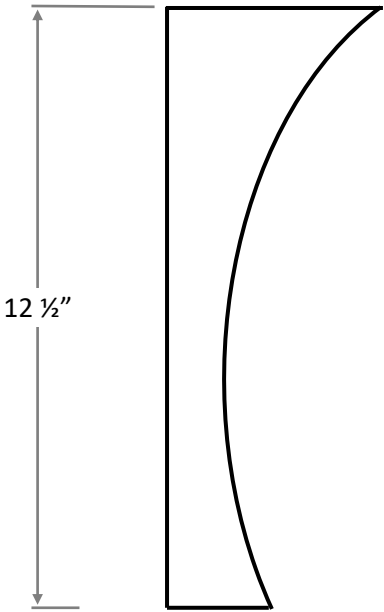



Outer Frame
4" = 1"

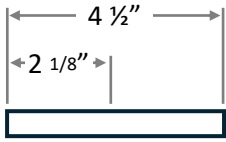




Attachment Arc for Load Cell Assembly
 $4'' = 1''$

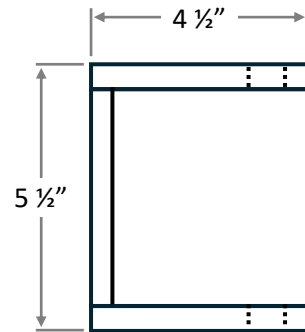


See detail for load cell assembly

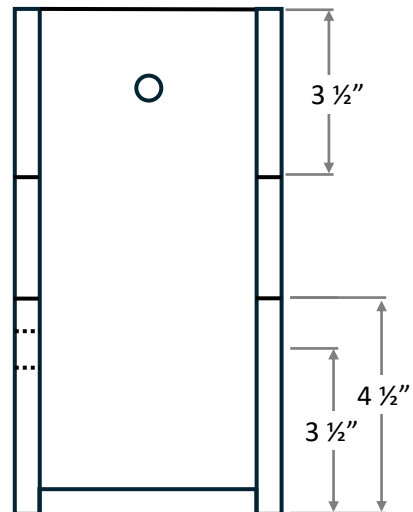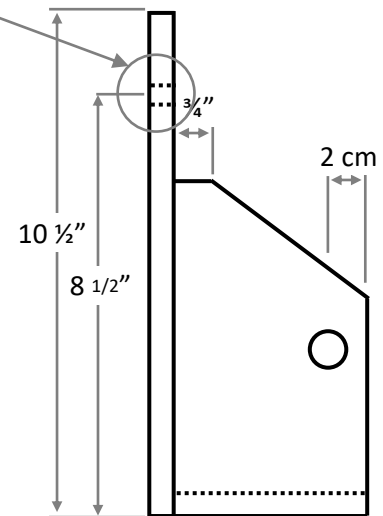



Foot Pedal
 $4" = 1"$

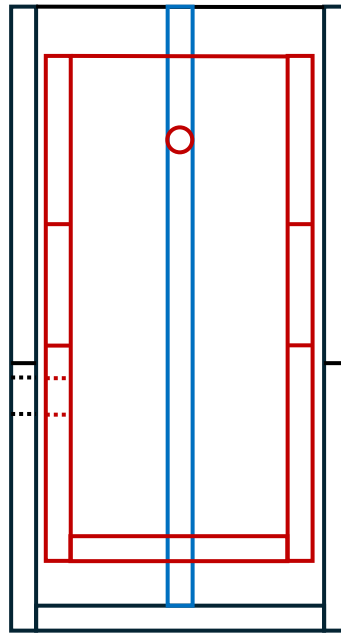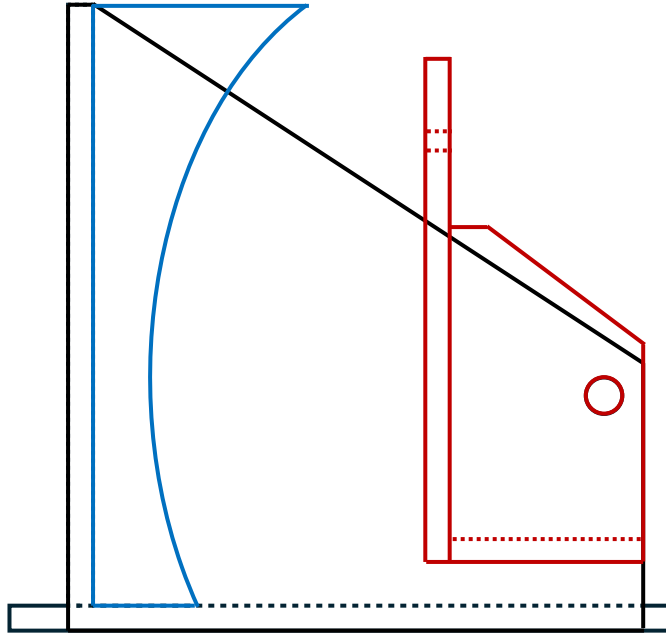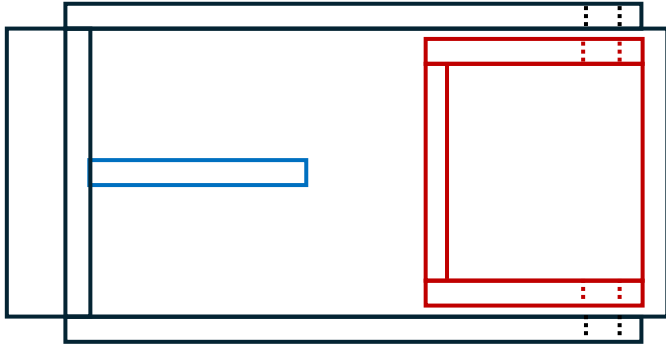

Entire Structural Assembly
 $4'' = 1''$

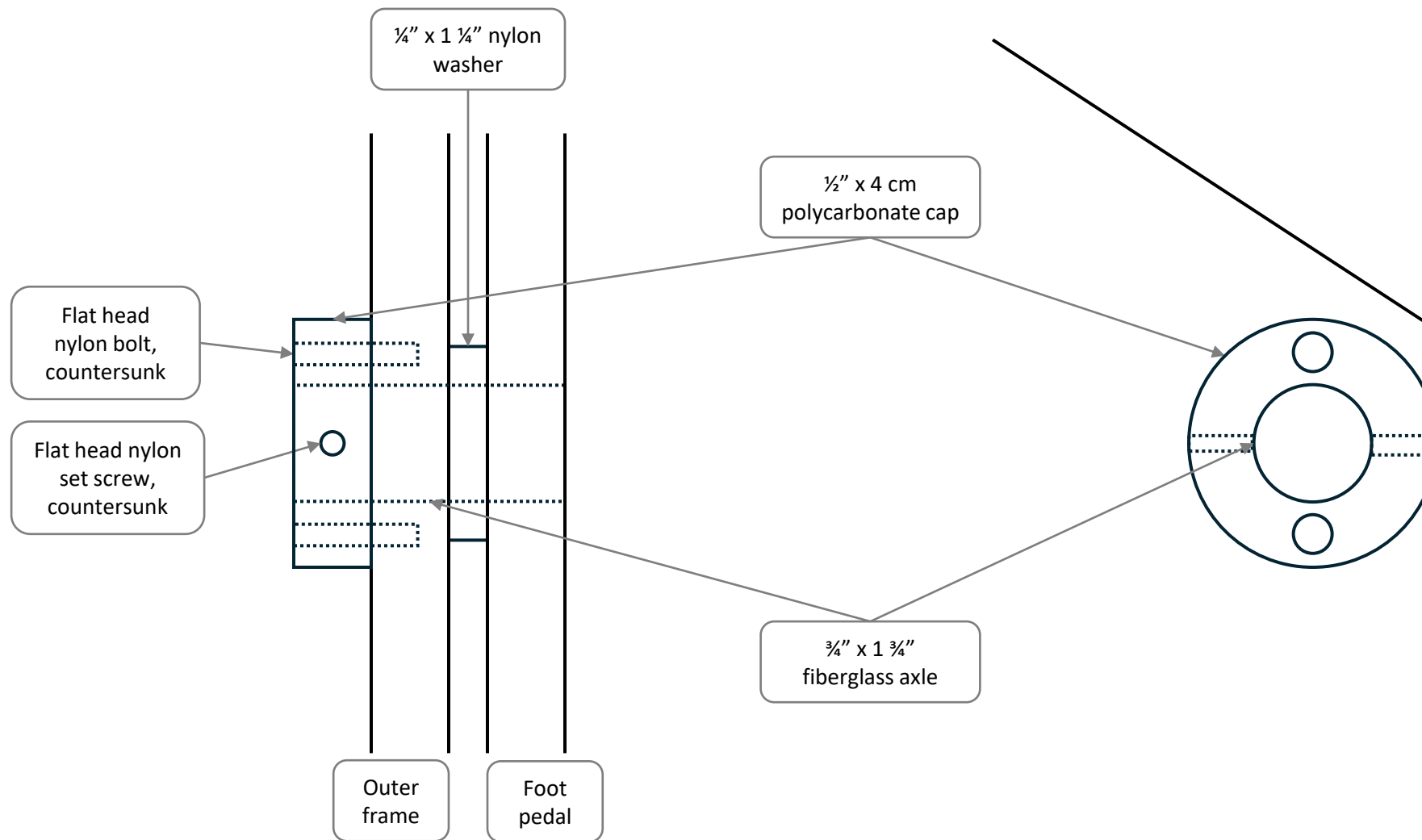


Axis Assembly
1" = 1"

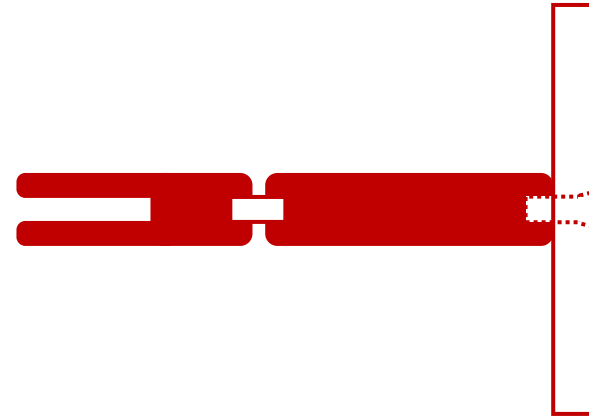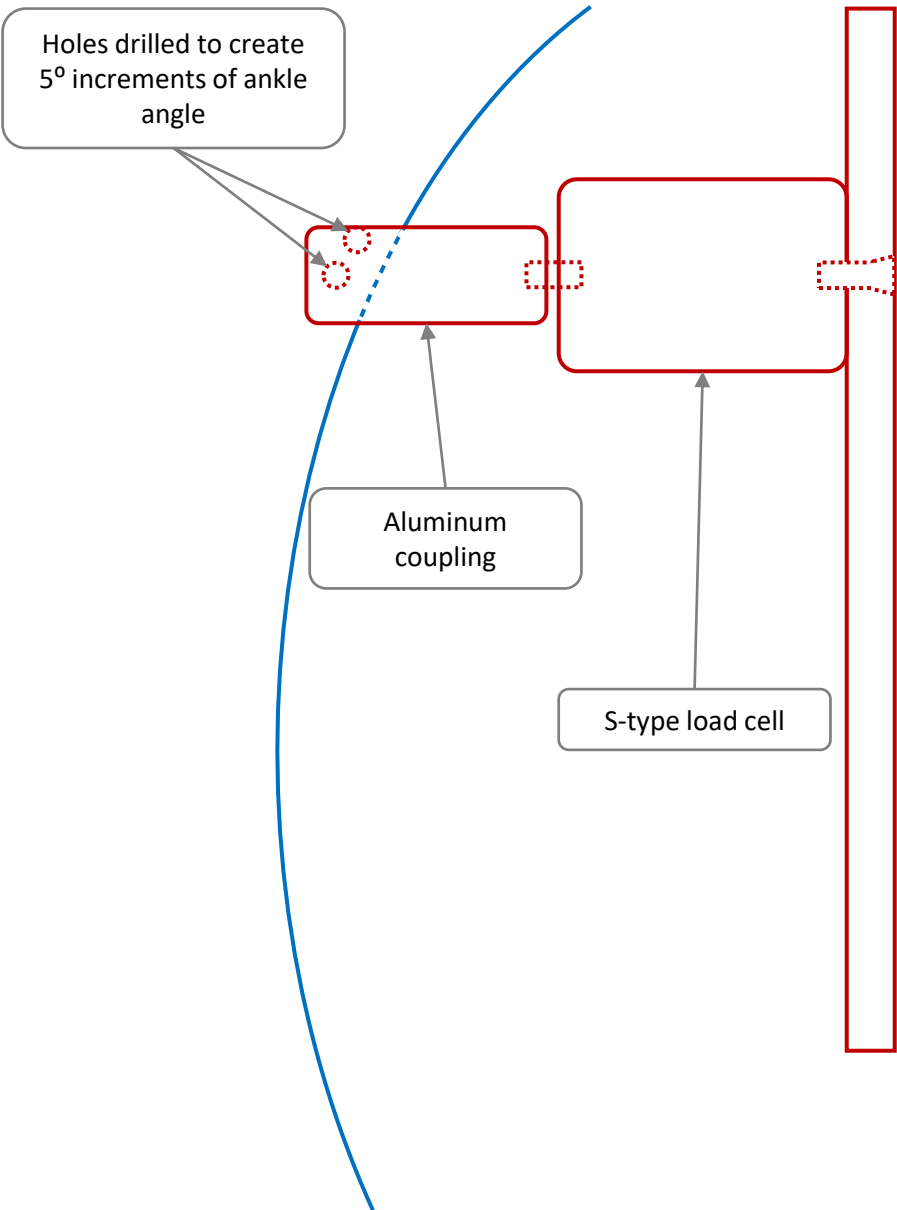

Verify radius and arc length of attachment arc, aluminum coupling size, bolt sizes and thread count, etc. as per load cell specifications

Secure coupling to attachment arc with cotter pin assembly


Load Cell Assembly
2" = 1"

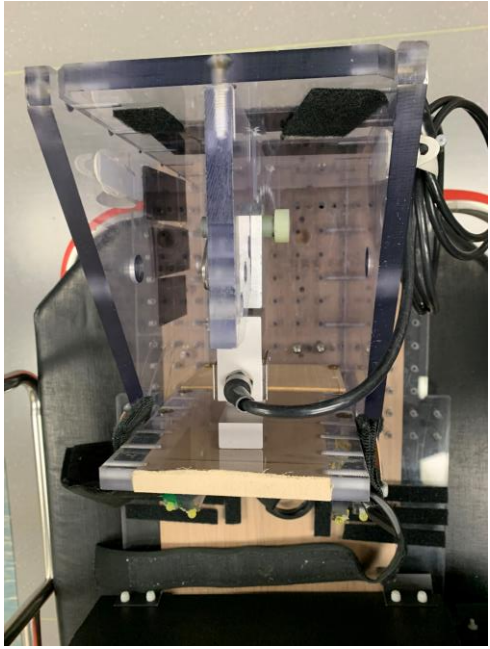

Top view




Oblique view

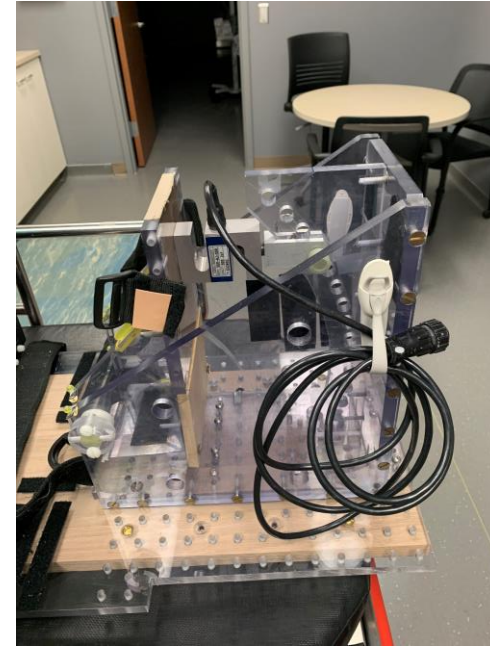

Side view




Load cell

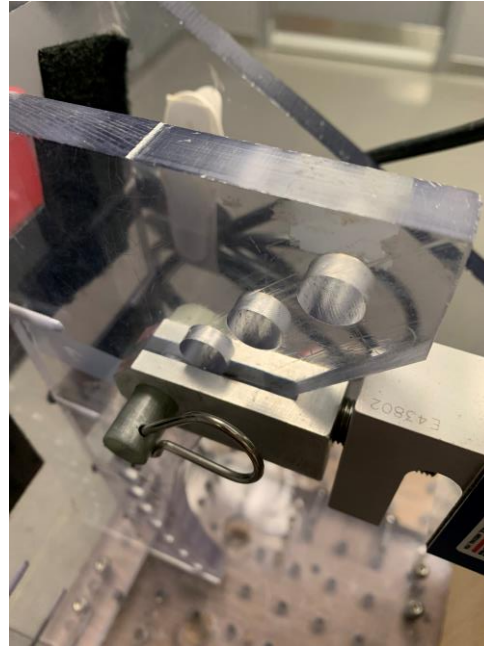

Load cell coupling to
frame



Electrical connection to
cable (cable passes
through the filter panel)

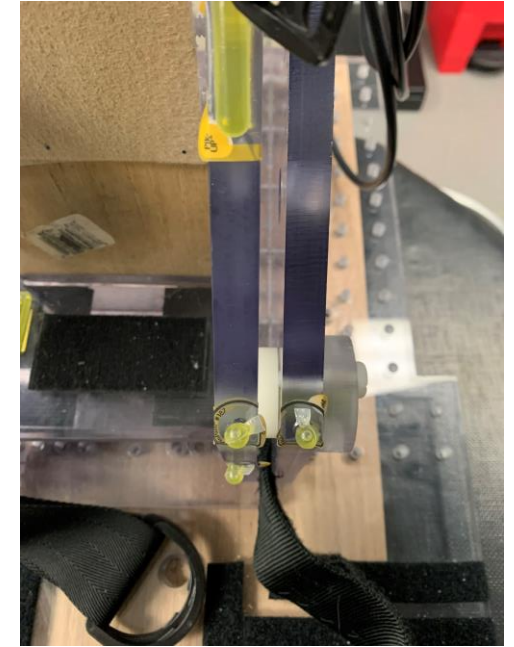

Axis assembly/rotary
joint for foot pedal
(Note fiducial markers for
alignment to
biomechanical
calculations)



Connection to force platform using stainless steel bolts



Long hex key for bolts



Strap over the ball of the foot



Strap over the ankle
(good for plantarflexion)